


How to Share Patient Data with Other Dentists




1

In DTX Studio™ Clinic:

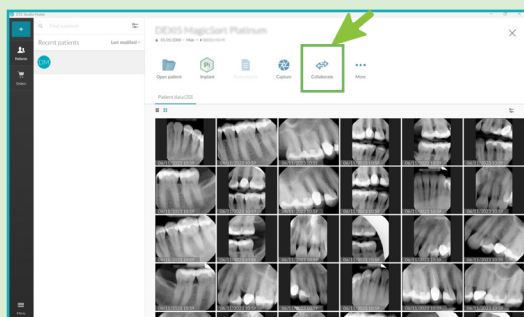


1. Select the patient case.
2. Click **Collaborate**.

2a

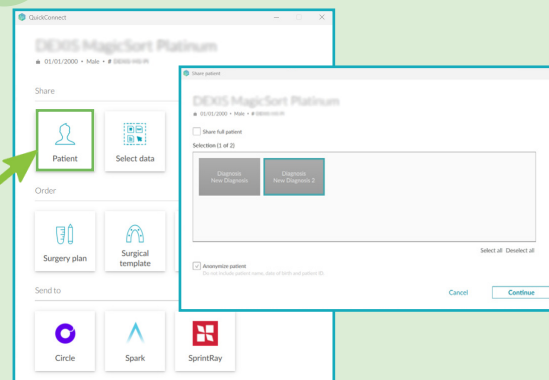

1. To share the whole case, click **Patient**.
2. Select the diagnoses you want to share, and click **Continue**.


2b

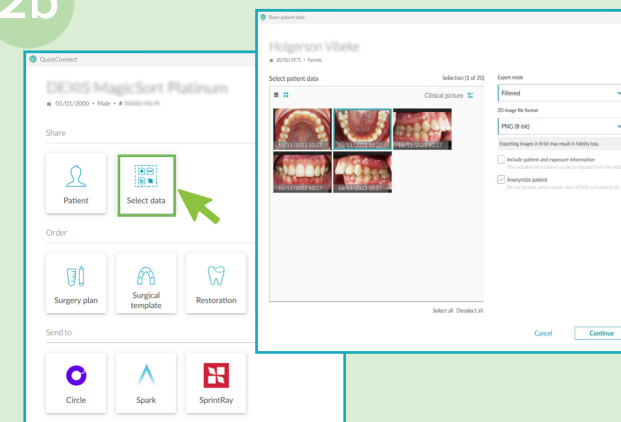

1. To share patient images, click **Select data**.
2. Select the images you want to share, and click **Continue**.


3

In DTX Studio™ Go:

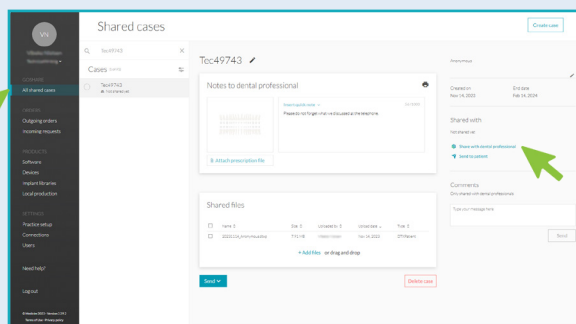


1. Complete the data in DTX Studio Go.
2. Click **Share with dental professional**.

4

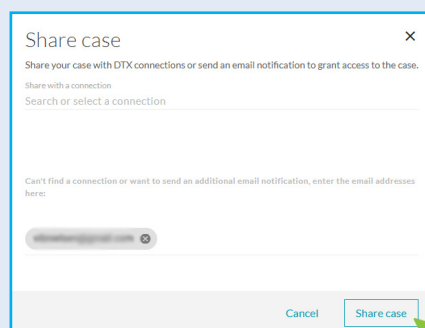

1. Enter the name of the connection or enter the email address.
2. Click **Share case**.

5

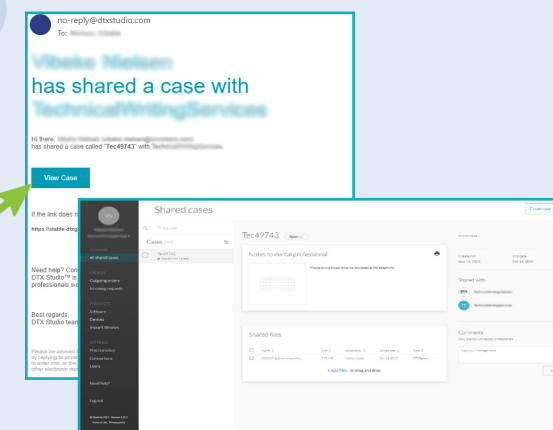

The partner receives an email with a link to the case on DTX Studio Go (**View case**).