



How to Create a Custom Report

1

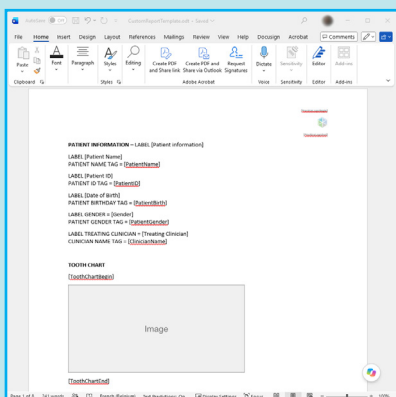


Design a template document in a word processor, such as Microsoft® Word, LibreOffice®, or Apache OpenOffice® Writer.

2

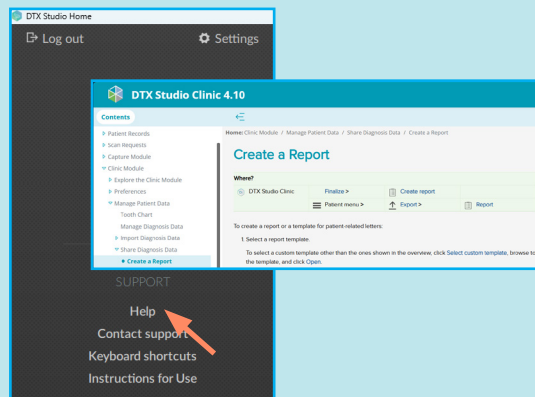


We recommend to start from an existing template. See [Report Template Example](#) for a sample template and 'Create a Report' in the help for available tags.

3

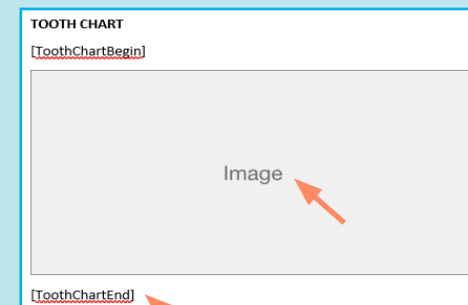


Close all tags correctly. Insert and resize placeholder images. These placeholders will be replaced with actual images when the report is created.

4

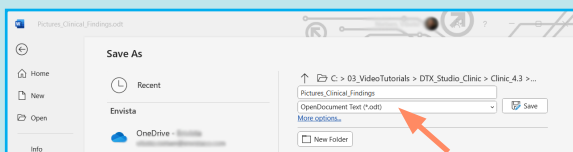

Save your template document as **OpenDocument Text (*.odt)**.


5

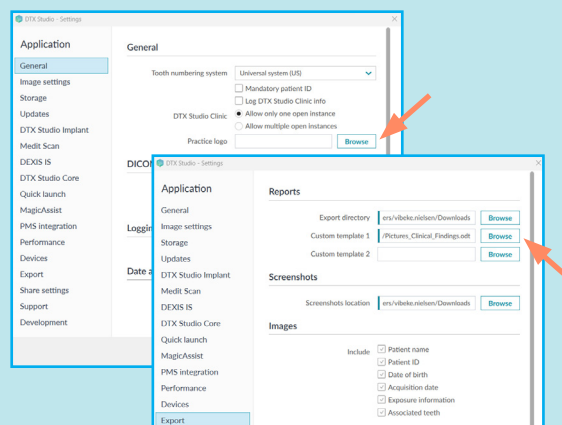

In the DTX Studio Clinic settings:

- Add your practice logo.
- Specify the location of up to two templates.

6

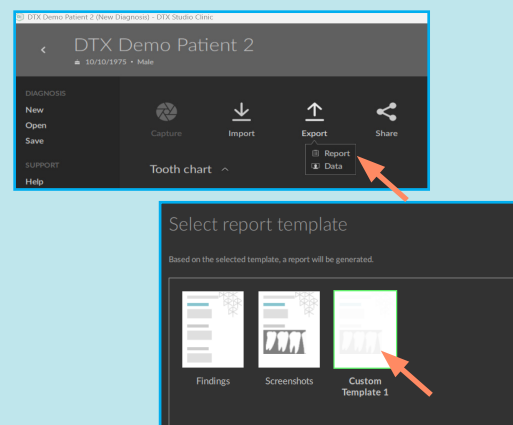

When generating a report for a patient in the Clinic module, select **Custom Template 1** (or 2).